

Chrom Tools For Everyone


Practical Tools for Better HPLC



Ferrule Puller Tool

Cat. No. 46025-FP


This nifty tool is very useful when a plastic, finger tight fitting or plug breaks off in a column or a female port. Stop throwing out columns or damaging them by removing the end fittings. 1 Each.

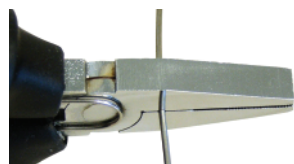

Beta Tool™ Tubing Straightner/Bender

Cat. No. 46010-BT

Bent tubing can cause internal kinks or simply turbulent flow. With this handy tool, you can remove kinks in tubing or make really smooth turns and bends in your SS tubing. 1 Each.

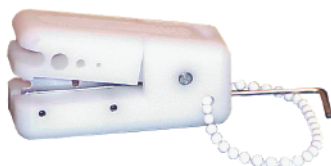

Tubing Cutter for PEEK

Cat. No. 49304-01

There are many low cost tubing cutters for PEEK tubing but if you want to easily make quality cuts in your tubing for best performance, you want to use this cutter. This tool has many features designed for quality cuts including a one sided blade to allow for 90° cuts and will hold your tubing so steady it is always a square cut. 1 Each.



The Valve Tool

Cat. No. 46090-V1

The Valve Tool has a standard 1/4" hex wrench on one side and a unique "open wrench" socket on the other. The open socket accommodates most standard male nuts with a hexagonal head. The "open section" of the tool allows you to place the tool on a male nut that from the top, over the tubing in the male nut which helps you tighten or loosen the installed nut. 1 Each.



HPLC Wrench

Cat. No. 46030-01

This wrench is 1/4" open end wrench on one side and 5/16" on the other side. The length is sufficient to get good leverage when you need to make sure your parts are tight enough. 1 Each.



Ferrule Removal Tool

Cat. No. 46020-FS

Use the Ferrule Removal Tool when you do not want to cut your stainless or Peek tubing to remove a damaged ferrule or to change the pilot depth when changing columns. 1 Each.

Other Tools are available, please inquire.

For more information visit www.mtc-usa.com

L-98



MTC • MicroSolv Technology Corporation

1 Industrial Way West, Bldg E, Eatontown, NJ 07724 USA

Ph. 1-732-578-1777

Fx. 1-732-578-9777

www.mtc-usa.com

Email: sales@mtc-usa.com