

Ketoconazole (1)

Channel B -- PDA Wavelength 214


345

345

Resolution (USP): 1.16397

| | |
|-------------|--|
| Buffer | CElixir pH 2.5 (25 gr) with HP- β -CD (2.5 gr) |
| Capillary | 50 μ x 40 cm |
| Separation | 25 kV |
| Injection | 5 seconds pressure |
| Temperature | 25°C |
| Instrument | Beckman P/ACE MDQ |

Au

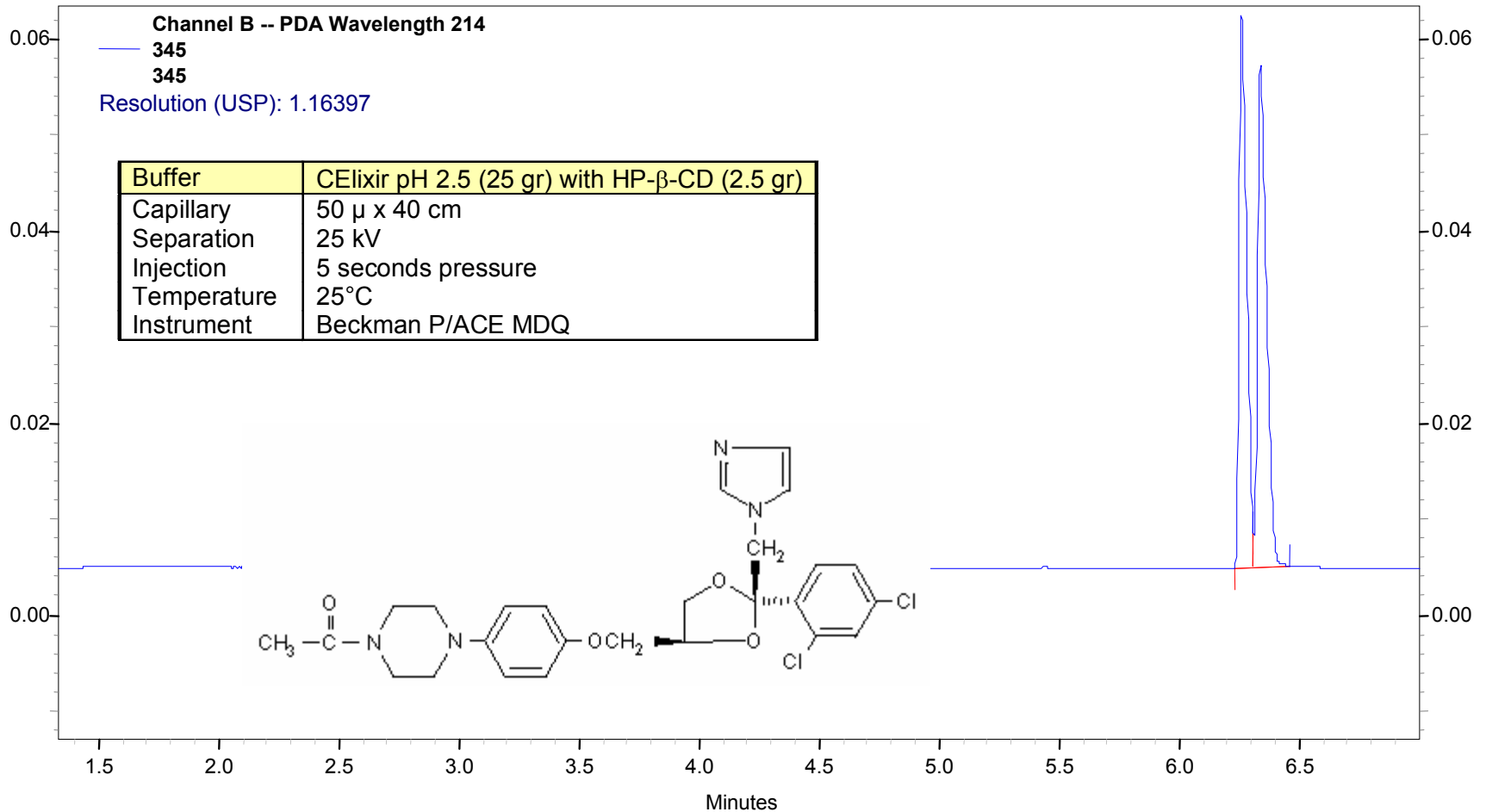

Au

Ketoconazole (1,2,3)

